For veterinary use only

## INTENDED USE

The BIONOTE Vcheck fPL 2.0 is an *in vitro* diagnostic test kit for the quantitative measurement of feline pancreas-specific lipase concentration in feline plasma or serum. This kit can be used for cats presenting with vomiting, fever, anorexia, depression or abdominal pain to assess the possibility of pancreatitis. The BIONOTE Vcheck fPL 2.0 is designed to be used only by veterinarians.

## PRINCIPLE

The Vcheck fPL 2.0 Test kit is a fluorescent immunoassay for the quantitative measurement of feline pancreas-specific lipase concentration. The Vcheck fPL 2.0 Test kit uses specific anti-fPL antibodies that bind to fPL. When the specimen is applied to the sample hole of the test device, fPL in the specimen migrates along the nitrocellulose membrane and forms complexes with the anti-fPL antibodies coated on the membrane. As a result, the density of the test line reflects the concentration of feline pancreas-specific lipase in the feline plasma or serum. The BIONOTE Vcheck Analyzer reads the density of this test line and calculates the feline pancreas-specific lipase concentration from the calibration curve data. The control line is a reference line which indicates the test has been performed correctly.

## MATERIALS PROVIDED

Reagent	5 Tests/Kit	10 Tests/Kit	20 Tests/Kit
① Vcheck fPL 2.0 Test device	5	10	20
② Assay diluent tube	5	10	20
③ Disposable pipette tip	10	20	40
④ Instructions for use	1	1	1

## MATERIALS REQUIRED, BUT NOT PROVIDED

- BIONOTE Vcheck Analyzer
- 25 µl pipette
- 100 µl pipette

## STORAGE AND STABILITY

- Store the test kit at 1~30 °C. **DO NOT FREEZE.**
- Do not store the test kit in the direct sunlight.
- The test kit is stable until the expiry date that is marked on the package label.

Reagent	Open status	Storage	Stability	Note
Test device	Unopened	1~30 °C, Sealed	15 months	Finished product
	Opened	Do not store	-	Use directly
Assay diluent	Unopened	1~30 °C, Sealed	15 months	Finished product
	Opened	Do not store	-	Use directly

## PRECAUTIONS

- The test kit is for feline use only. Do not use for other animals.
- The test device is sensitive to humidity and heat. Perform the test immediately after removing the test device from the aluminum foil pouch.
- Do not reuse test components.
- Do not touch the membrane in the result window of the test device.
- Do not use the test kit beyond the stated expiry date marked on the label.
- Do not use the test kit if the pouch is damaged or the seal is broken.
- Do not mix components from different lot numbers, the components in this kit have been quality control tested as a standard batch unit.
- All samples should be handled as being potentially infectious. Wear protective gloves while handling samples. Wash hands thoroughly afterwards.
- Decontaminate and dispose of all samples, used kits and potentially contaminated materials safely in accordance with national and local regulations.

- BIONOTE Vcheck analyzer is recommended to use at 15~30 °C.
- Severely hemolyzed samples or serum with precipitate such as fibrin may give erroneous results.
- A professional veterinarian should make a final diagnosis based on the results of this product, other test results and clinical findings.
- Strictly follow the test procedure (E.g. adequate sample volume), as failure to do so may adversely affect test performance and/or produce invalid results.
- Although the Vcheck fPL Test Kit offers simple and quick quantitative measurement of fPL concentration in feline serum or plasma, there may be a difference in the detection performance with other clinical or laboratory methods with more sophisticated principles.

## COLLECTION AND PREPARATION OF SAMPLE


- Feline plasma or serum should be used with this test.
- [Plasma]**
  - Blood should be collected with a disposable syringe and added to a tube containing **anticoagulant (EDTA)**, and then separate plasma by centrifugation.
  - If plasma samples are not tested immediately, they should be refrigerated at 2~8 °C and used within 2 weeks. For longer storage, plasma can be frozen (-20 °C or colder). Frozen samples should be brought to room temperature (15~30 °C) prior to use.
- [Serum]**
  - Collect the whole blood into a blood collection tube containing NO anticoagulant. Leave to settle for 30 minutes for blood coagulation and then centrifuge to obtain a serum supernatant.
  - If serum samples are not tested immediately, they should be refrigerated at 2~8 °C and used within 2 weeks. For longer storage, serum can be frozen (-20 °C or colder). Frozen samples should be brought to room temperature (15~30 °C) prior to use.

## TEST PROCEDURE

All reagents and samples must be at room temperature (15~30 °C) before use.

### [Coding]

- Turn on V200 Analyzer and select "Standard Test".

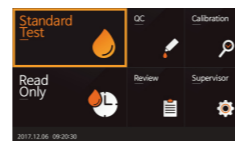
- Remove the test device from the aluminum foil pouch. Once the "Insert Device" is displayed in the screen, insert the test device.


### [Dilution of sample & Measurement]


- Use a 25 µl pipette to draw of serum (plasma) 25 µl and add the sample into an assay diluent tube (100 µl).

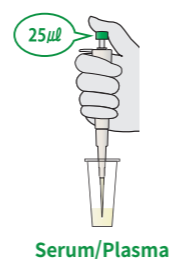- Use a 100 µl pipette to mix the sample with diluent by pipetting for 5~6 times.

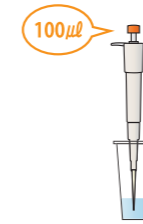
- Add the mixed sample (100 µl) into the sample hole of the test device. Then, press the "Start" button.

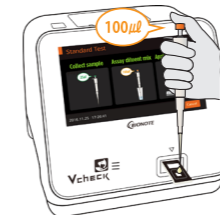- If the time to press [START] button is delayed, it may affect the test result.

- The V200 analyzer will display the test result on the screen after 15 minutes.

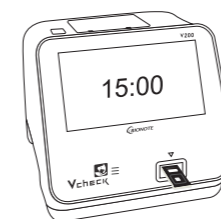- Remove the test device.

\* **Strictly follow the test procedure including the amount of sample (25 µl) used and the test time (5 min), as failure to do so may adversely affect test performance and/or produce invalid results.**

## INTERPRETATION OF THE RESULT

- Read the concentration value of feline pancreas-specific lipase appearing on the display of the BIONOTE Vcheck Analyzer. (1~50 ng/mL)
- If "↓ 1 ng/mL" appears on the display, it means the concentration of feline pancreas-specific lipase in the specimen is less than 1 ng/mL.
- If "↑ 50 ng/mL" appears on the display, it means the concentration of feline pancreas-specific lipase in the specimen is greater than 50 ng/mL.
- If the [Invalid] result appears on the screen, a retest shall be carried out.

## REFERENCE RANGE

≤ 3.5 ng/mL	Normal
3.6 ~ 5.3 ng/mL	Suspected
≥ 5.4 ng/mL	Consistent with pancreatitis

\* It is recommended that each laboratory establishes its own reference range.

## PERFORMANCE CHARACTERISTIC

- Measuring Range**  
fPL concentration can be measured within the range of 1~50 ng/mL by using feline serum or plasma. To get numerical results for ↑ 50 ng/mL concentration ranges, samples need to be diluted with negative serum and measured again. To calculate the final fPL concentration, the measured result must be multiplied by the dilution factor.
- Interference**  
No interference was observed for each substance up to the concentration presented in the following table.

Interfering substances	Concentration
Hemoglobin	< 10,000 mg/dL
Bilirubin	< 5 mg/dL
Triglyceride	< 10 mg/dL
Cholesterol	< 10 mg/dL
Vitamin C	< 20 mg/dL

## SCREEN MESSAGES AND TROUBLE SHOOTING

[V200]

Error message	Error description
Contaminated Device	The test device is damaged or inserted improperly. <b>Solution:</b> Discard the test device and retest with a new test device and a new specimen.
Insufficient Sample	An insufficient amount of blood has been applied. <b>Solution:</b> Retest with a new test device with enough specimen, ensuring that blood is placed in to the narrow channel in the top edge of the test device.
Expired Device	The test devices are expired. <b>Solution:</b> Retest with a new test device that is not expired.
Temperature Error	The environmental temperature is above or below the operating range of the analyzer. <b>Solution:</b> Move to an area in the acceptable temperature range for the analyzer and perform the test. Do not heat or cool the analyzer artificially.
Printer Connection Fail	The communication between analyzer and barcode or printer has failed. <b>Solution:</b> Reconnect the analyzer and external device. If the error continues after turning ON/OFF the analyzer, please contact BioNote, Inc.
Barcode Error	The measured total hemoglobin is out of the range of 7 to 23g/dL. <b>Solution:</b> This error occurs when a specimen has a total hemoglobin in the abnormal range. If the error continues after turning ON/OFF the analyzer, please contact BioNote, Inc.
Extremely High Total Hemoglobin	The test is invalid. <b>Solution:</b> Retest with a new test device and a new patient specimen. If the error continues after turning ON/OFF the analyzer, please contact BioNote, Inc.
Result: Invalid	The calibration is overdue. <b>Solution:</b> If the error continues after turning ON/OFF the analyzer, please contact BioNote, Inc.
Calibration Overdue	Loading a test device that is not supported by the analyzer. <b>Solution:</b> Check whether the test device is manufactured by BioNote, Inc.
Not Supported Device	Internal error has occurred. <b>Solution:</b> If the error continues after turning ON/OFF the analyzer, please contact BioNote, Inc.
EEE	

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